


# BetterLiving

## Pressure Relief Cushion

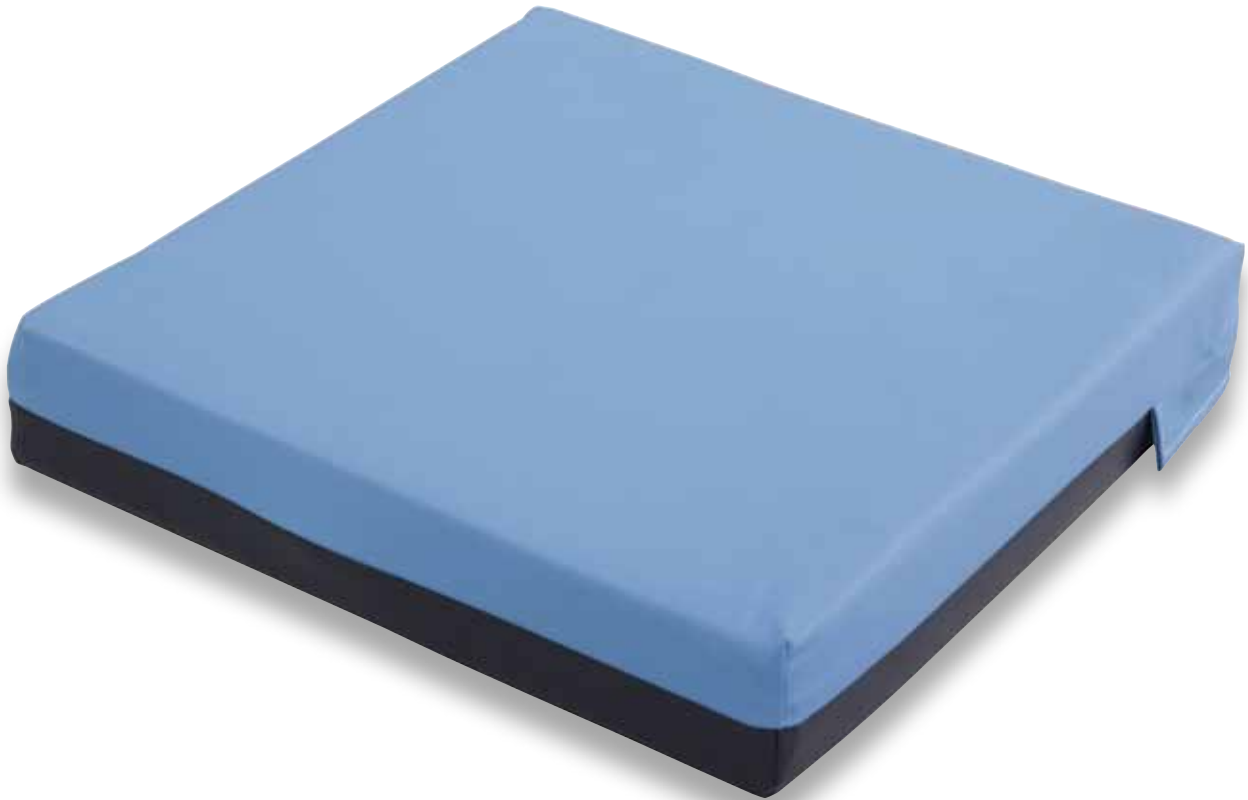

This two layer foam cushion provides improved pressure redistribution when seated.

A two layer foam cushion for improved pressure redistribution when seated. This multi-purpose product is suitable for users with mild postural needs and comprises a high density foam support base with a super soft memory foam top layer. Enclosed by a zippered two-way stretch waterproof polyurethane fabric cover for enhanced anti-shear and anti-friction benefits, with 360° zip and waterfall skirt to protect from moisture ingress. Seams are stitched for strength, then welded for or maximum infection control.

› Formerly called the Comfort Cushion

Item	Pressure Relief Cushion
Code	WSGS7258
Height	75 mm
Width	410 mm
Depth	410 mm
Capacity	100 kg
Warranty	1 year